



P/V Tool® Pro

Recruitability assessment and lung recruitment

The Protective Ventilation Tool (P/V Tool Pro) performs a respiratory mechanics maneuver that records a quasi-static pressure/volume curve. This method can be used to assess lung recruitability and determine the recruitment strategy to apply.

The P/V Tool Pro can also be used to perform a sustained inflation recruitment maneuver and measure the increase in lung volume. It is particularly helpful for ARDS patients, as selecting an appropriate lung-recruitment strategy and the correct PEEP setting as an anti-derecruiting force can be critical for this patient group.

- ✓ Tool for individualized, lung-protective ventilation
- ✓ Simple and safe way to assess patients' potential for recruitment and perform lung recruitment maneuvers
- ✓ No need to disconnect the breathing circuit or make changes to the mode or ventilator settings
- ✓ Easily repeatable process to monitor changes in patient condition and effectiveness of treatment over time
- ✓ Interpretation is aided by automatic calculations and cursors to assist with analysis

The P/V Tool Pro is available as an option on the HAMILTON-G5, HAMILTON-C6, and HAMILTON-C3 mechanical ventilators and is a standard feature on the HAMILTON-S1.



Swiss
Quality

HAMILTON
MEDICAL
Intelligent Ventilation since 1983

Assess recruitability and perform a recruitment maneuver

A diagnostic tool: Assess recruitability of ARDS patients (graphic 1)

Using the pressure/volume curve, it is possible to distinguish patients with low recruitability, in whom recruitment maneuvers and high PEEP are not appropriate, from patients with high recruitability, who may benefit from recruitment maneuvers and higher PEEP.¹


Assessing the lung recruitability is a prerequisite for determining the optimal recruitment strategy for an ARDS patient.²

A recruitment tool: Perform a recruitment maneuver (graphic 2)


The P/V Tool Pro helps you perform a safe, sustained inflation recruitment maneuver. The pressure ramp, maximum pressure, duration, and final pressure (PEEP level after the recruitment maneuver) can be individualized for each patient.

The P/V Tool Pro offers you a fully controllable and repeatable method to perform a recruitment maneuver, and keeps track of how effective the maneuver is.

1 Demory D. Intensive Care Med 2008; 34:2019–2025. | 2 Gattinoni L. N Engl J Med 2006;354:1775-86.

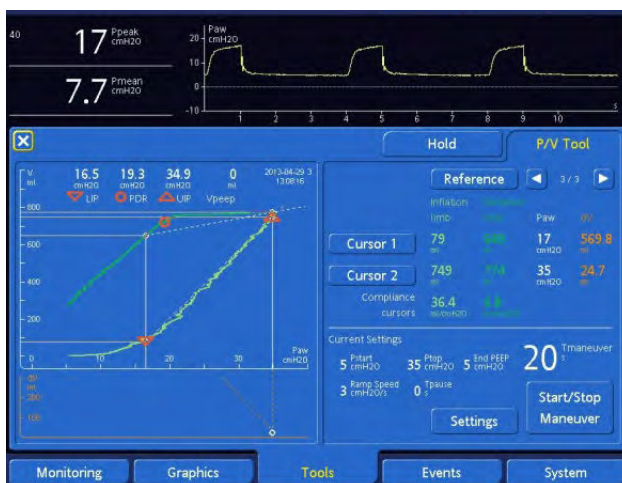

Graphic 1: A cursor function enables graphical analysis of the curve, including identification of inflection points

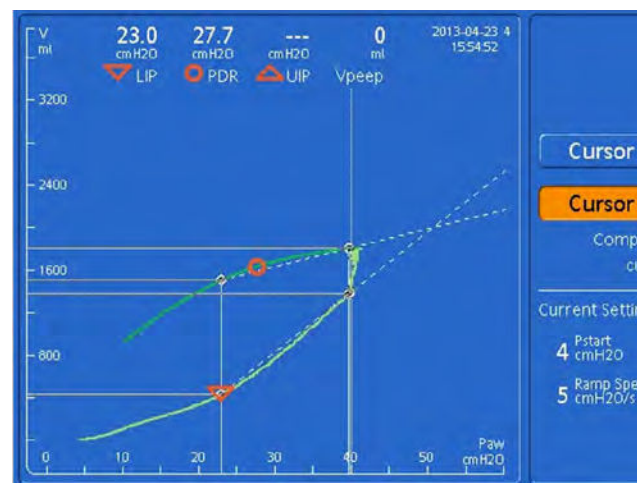

Graphic 2: The increase in lung volume resulting from the recruitment maneuver can easily be measured using the cursors

Apply lung-protective strategies

A lung-protection tool

Lung protection is a major objective for every patient on mechanical ventilation. Lung collapse induces regional inhomogeneity and increases the risk of atelectrauma.³ A well-conducted recruitment strategy, combining recruitment maneuvers and adequate PEEP, increases lung homogeneity.

Lung protection not only reduces mortality in ARDS patients, but also reduces the risk of secondary ARDS in normal lung patients, as well as the complications in post-surgery patients.^{4, 5, 6,}

An advanced tool combined with esophageal pressure monitoring

In combination with esophageal pressure measurement, the P/V Tool Pro can give you a clearer understanding of the lung and chest-wall mechanics. This enables you to apply a lung-protective ventilation strategy by titrating the PEEP level⁷, and optimizing the parameters for the recruitment maneuver, driving pressure and tidal volume.

3 Caironi P. Am J Respir Crit Care Med 2010; 181: 578–586. | 4 Acute Respiratory Distress Syndrome Network. N Engl J Med. 2000 May 4;342(18):1301-8. | 5 Determann RM. Crit Care. 2010;14(1):R1. doi: 10.1186/cc8230. Epub 2010 Jan 7. | 6 Futier E. N Engl J Med 2013;369:428-437. | 7 Talmor D. N Engl J Med. 2008 Nov 13;359(20):2095-104. doi: 10.1056/NEJMoa0708638. Epub 2008 Nov 11.

“

We use the P/V Tool Pro to optimize initial PEEP settings on almost all ventilator patients once they have been sedated for the intubation. We also use the recruitment portion of the P/V Tool Pro a lot, particularly for patients who are having repeated atelectasis.

Ken Hargett, Director Respiratory Care
Houston Methodist Hospital, Houston (TX), USA





“

We recommend that our respiratory therapists use the P/V Tool Pro to optimize PEEP as soon as they put the patient on the ventilator. The therapists find it very helpful, especially on the sicker patients.

Camille Neville, Critical Care Educator Respiratory Care
Florida Hospital, Orlando (FL), USA

Manufacturer:

Hamilton Medical AG

Via Crusch 8, 7402 Bonaduz, Switzerland

☎ +41 (0)58 610 10 20

info@hamilton-medical.com

www.hamilton-medical.com

689202.04

Specifications are subject to change without notice. Some features are options. Not all features are available in all markets. For all proprietary trademarks (®) and third-party trademarks (®) used by Hamilton Medical AG see www.hamilton-medical.com/trademarks. © 2018 Hamilton Medical AG. All rights reserved.

P/V Tool Pro