

# HAMILTON-HF90

## ROX index



### What is the ROX index?

When delivering high flow oxygen therapy (HFOT), recognizing when to escalate to noninvasive ventilation (NIV) or endotracheal intubation (EIT) is critical. The ROX index (Respiratory rate rate – OXYgenation) is an easy-to-use tool that can aid clinicians in making this decision at pivotal moments during HFOT. The index relies on parameters directly linked to oxygenation (assessed by SpO<sub>2</sub> and FiO<sub>2</sub>) and respiratory distress (assessed by RR). The respiratory rate can be monitored using Masimo RRp<sup>1</sup> technology.


### How to calculate and interpret it

The ROX index is calculated by dividing the oxygen saturation (SpO<sub>2</sub>) by the fraction of inspired oxygen (FiO<sub>2</sub>) and then dividing the result by the respiratory rate (RR):  $\text{ROX index} = (\text{SpO}_2 / \text{FiO}_2) / \text{RR}$ .

To measure the ROX index, values for SpO<sub>2</sub>, FiO<sub>2</sub>, and RR are recorded and inserted into the formula. A higher ROX index indicates better respiratory function and effectiveness of high flow oxygen therapy, while a lower index signals a potential deterioration in respiratory function. However, regular monitoring and trend analysis are also crucial, as they are more informative than a single measurement.<sup>1</sup>

ROX index calculation	Example I		Example II		Result description
ROX index = (SpO <sub>2</sub> / FiO <sub>2</sub> ) / RR	ROX ≥ 4.88 success		ROX < 3.85 consider intubation		ROX index ≥ 4.88 after two hours of treatment, indicates a high probability that intubation will not be necessary.
	SpO <sub>2</sub>	94%	SpO <sub>2</sub>	92%	
	FiO <sub>2</sub>	0.6	FiO <sub>2</sub>	0.8	
	SpO <sub>2</sub> / FiO <sub>2</sub>	157	SpO <sub>2</sub> / FiO <sub>2</sub>	115	ROX index < 3.85 indicates a higher risk of treatment failure.
	Rate	25	Rate	35	
	ROX	6.27	ROX	3.29	

## Recommendations for measuring the ROX index

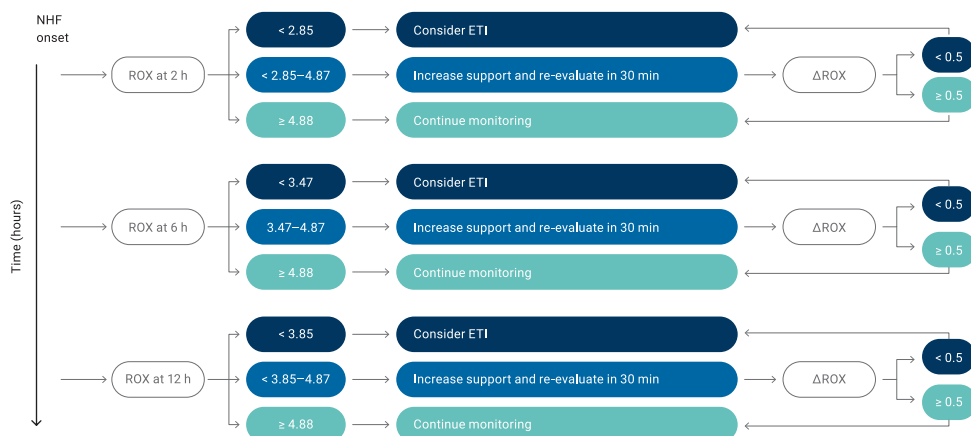


## Accessing ROX monitoring on the HAMILTON-HF90

Select Trends

Tap the right arrow

Tap until ROX index appears

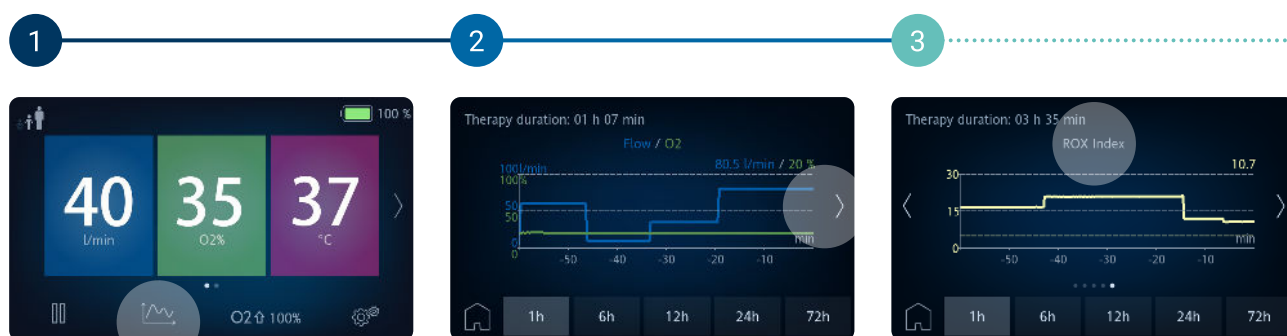

## Comparative insights in non-COVID-19 and COVID-19 pneumonia

In acute hypoxemic respiratory failure due to non-COVID-19 pneumonia, the ROX index identified patients at low risk of HFOT failure with a cut-off value of 4.88 measured after 12 hours of HFOT.<sup>1</sup>

In COVID-19 patients, five retrospective studies<sup>2,3,4,5,6</sup> showed that patients with a successful outcome had a higher ROX index, but the cut-off point for values associated with success varied between 5.55 after 6 hours<sup>4</sup> and 3.67 after 12 hours<sup>5</sup>.



Hamilton Medical AG

Via Crusch 8, 7402 Bonaduz, Switzerland

+41 (0)58 610 10 20

info@hamilton-medical.com

www.hamilton-medical.com



Scan the QR code for the list of references.

The products shown here are not available for purchase by the general public. The information provided here is intended for healthcare professionals only. Always read the labels and follow the product's instructions for use. Specifications are subject to change without notice. Some features are options. Not all features are available in all markets. All images are used for illustrative purposes only and may not accurately represent the product or its use. For all proprietary trademarks (®) and third-party trademarks used by Hamilton Medical AG see [www.hamilton-medical.com/trademarks](https://www.hamilton-medical.com/trademarks). ©2024 Hamilton Medical AG. All rights reserved.