



**Laerdal®**  
helping save lives



# Bystander CPR Saves Lives

## Little Anne® QCPR

Providing feedback based on just watching students perform can be challenging for instructors. The next generation of our Little Anne® manikin includes QCPR measurement and feedback technology that takes the guesswork out of instruction—all from the palm of your hand. Little Anne QCPR works with an Instructor app that allows you to monitor and coach up to six students at once, a Learner app that lets students monitor their own performance, and a Classroom app for training large groups.

Little Anne QCPR helps instructors improve CPR training quality and efficiency in their courses, while adding a fun and engaging element for students that promotes learning and motivation¹.



**laerdal.com**



# Little Anne® QCPR

Little Anne® QCPR is still the same durable, lightweight manikin as the original Little Anne, but now includes QCPR technology that provides:

- Real-time feedback on compressions and ventilations
- A summary with an overall score and improvement suggestion
- Details on compression release, depth and rate, ventilation volume, number of compressions/ventilations/cycles

## Feedback Options

When used with your choice of feedback option(s), Little Anne QCPR meets the AHA's requirement for feedback devices.



### QCPR Instructor app

*For instructors to monitor multiple students*



- Free download on App Store and Google Play
- Connects to 1-6 manikin(s)
- Shows live and summative feedback of CPR performance
- Provides suggestion for improvement, post-session
- Results are saved in the app
- QCPR Race game reinforces quality CPR skills in a fun and engaging way




### QCPR Learner app

*For students to monitor their own performance*



- Free download on App Store and Google Play
- Connects to 1 manikin
- Shows live and summative feedback of CPR performance
- Provides suggestion for improvement, post-session
- Results are saved in the app

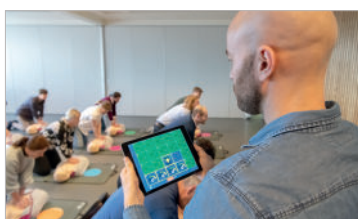

### QCPR Classroom app

*For instructors to monitor large groups*



- Free download on App Store - iPad only
- Connects to 1-42 manikin(s)
- Shows live and summative feedback of CPR performance
- Provides group performance summary
- QCPR Race game reinforces quality CPR skills in a fun and engaging way



### SkillGuide™

- Provides a tethered, "plug and play" option
- Can be used standalone or as a supplement to either of the apps

Little Annes manufactured after February 1999 can be quickly and easily upgraded to QCPR with the Little Anne QCPR Upgrade Kit.



## Ordering Information

123-01050	Little Anne QCPR
123-03050	Little Anne QCPR Dark Skin
124-01050	Little Anne QCPR 4-pack
124-03050	Little Anne QCPR 4-pack Dark Skin
123-60750	Little Anne QCPR Upgrade Kit – compatible with Little Annes manufactured after Feb 1999

## Accessories and Consumables

123-30050	SkillGuide with extension cable
020300	Little Anne Airways (pkg. 24)
020301	Little Anne Airways (pkg. 96)
15120103	Manikin Face Shields (6 rolls of 36 ea)
152400	Manikin Wipes (pkg. 50)
152401	Manikin Wipes (pkg. 1200)
310210	Manikin Faces (pkg. 6)
310215	Manikin Faces, Dark Skin (pkg. 6)
123-40050	QCPR Race Circles, 6-Pack

Visit [Laerdal.com/LittleAnneQCPR](http://Laerdal.com/LittleAnneQCPR) to learn more about Little Anne QCPR and how you can upgrade your existing Little Anne manikins.

©2018 Laerdal Medical. QCPR is a registered trademark of Laerdal Medical. All rights reserved. Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries. Google Play and the Google Play logo are trademarks of Google Inc. #18-17163

I. Smart et al. Does real-time objective feedback and competition improve performance and quality in manikin CPR training – a prospective observational study from several European EMS. Scandinavian Journal of Trauma, Resuscitation and Emergency Medicine. 2015



**Laerdal®**  
helping save lives