

# microlife®

A partner for people. For life.

Compact and  
**PORTABLE**

Small, silent and  
light solution for the  
whole family



## NEB800 Mesh Nebuliser

- Portable, pocket size design
- One button – User friendly interface
- Ultra light and silent
- Nebulisation rate: 0.2 ml/min. (NaCl 0.9%)
- Particle size:  $\approx 5 \mu\text{m}$  (MMAD)
- Acoustic noise level: < 30 dBA

### Key features




Portable – Easy for travel



Ultra light and silent



Biocompatible  
Phthalate and BPA free material

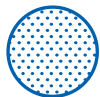

Stainless steel mesh




Auto off

Automatic switch-off

### Efficient Nebulisation System

- 1 Medication (prescription only).
- 2 Aerosol of medication:  $5 \mu\text{m}$ , capable of penetrating the drug quickly and deeply into the lungs.
- 3 Efficient venturi system: Draws in ambient air, quickly delivers the medication particles, so it shortens the therapeutic session.
- 4 Mesh:
  - High anti-oxidized mesh delivers stable particles whilst providing a longer lifetime for the device.
  - High vibrating frequency (approximately 103 ~ 123 kHz) increases the nebulisation rate.

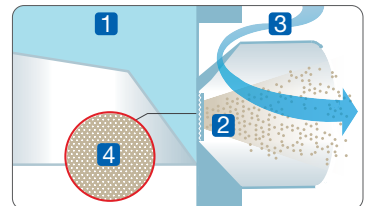

### Fully equipped

- With adult and child face masks, mouthpiece and carrying bag
- Easy storage and transportation
- Intuitive design for comfortable filling of medication



CE 1639

2 year  
guarantee