



MADgic Airway

Intubating Airway with Mucosal Atomization and Oxygen Delivery

Difficult and Awake

The MADgic Airway Device from Teleflex combines atomized topical anesthetic and oxygen delivery in an innovative and elegantly designed fiberoptic compatible oral airway.

Exceptional Coverage

The MADgic Airway Device protects the fiberoptic scope, and allows passive oxygenation with the option to place mask over the patient to deliver additional oxygen or ventilation. The built-in atomizer easily reaches the vocal cords and larynx and trachea with a consistent, fine atomized spray.

Versatile Design

The device conforms to patient anatomy allowing topical anesthetic to be delivered to the patient while in an upright or recumbent position, and while awake or asleep. The MADgic Airway Device provides a single size solution that allows for a wide range of endotracheal tubes, giving the ability to easily side unload ETTs and be removed after ETT insertion.

Compatible

The MADgic Airway Device may be used in conjunction with the MADgic Laryngo-Tracheal Mucosal Atomization Device* for difficult and awake intubation and extubation cases requiring a fiberoptic scope.

Benefits



The Clinician

Provides an innovative solution to create a temporary airway to facilitate fiberoptic intubation



Your Institution


Allows for preoxygenation, passive oxygenation, and topical anesthetic application with the cost of a single intubating airway



The Patient

Increase patient comfort and reduced risk of patient injury

*For use with drugs approved for intranasal and oropharyngeal delivery.

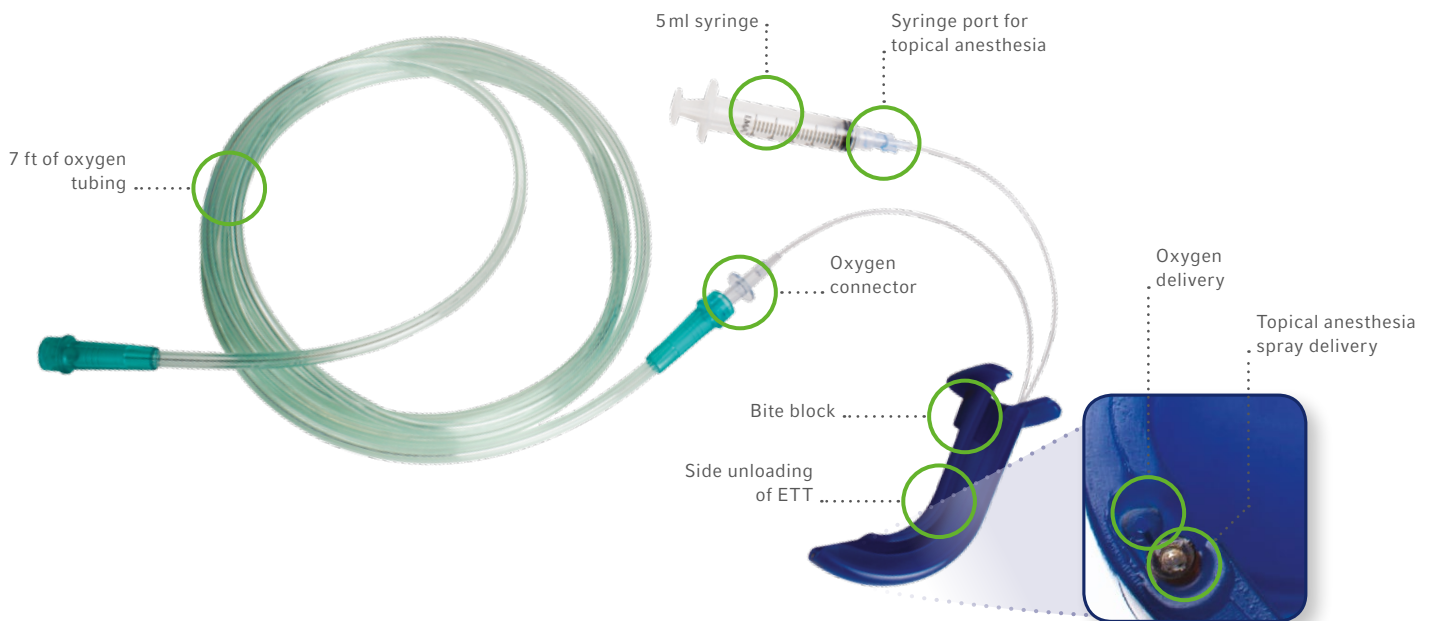

MADgic Airway Device Specifications

Typical Particle Size	30–100 microns
System dead space	0.15 ml
Oxygen Flow Rate	2–3 l/min @ 50 psi
Size	9 cm Airway (6.5–8.0 ETT)

MADgic Airway Device

ITEM NUMBER	PRODUCT DESCRIPTION	QTY/BOX
MAD780	MADgic Airway Device with 5 ml syringe and 7 ft long green oxygen tubing	10

The MADgicWand Pharyngeal Tissue Retractor/Atomizator Accessory

The MADgicWand Pharyngeal Tissue Retractor/Atomizer from Teleflex allows thorough topicalization and precise anesthetic placement to the difficult-to-reach glottic opening.



MADgicWand Device Specifications

Typical Particle Size	30–100 µm
System dead space	0.24 ml

MADgicWand Pharyngeal Tissue Retractor/Atomizer

ITEM NUMBER	PRODUCT DESCRIPTION	QTY/BOX
MAD650	MADgicWand Device with 5 ml Syringe	20

To learn more, visit Imaco.com/products

Teleflex, the Teleflex logo, MAD, MADgic, MADgic Airway and MADgicWand trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this document is not a substitute for the product Instructions for Use. The products in this catalogue may not be available in all countries. Please contact your local representative. All data current at time of printing (07/2016). Subject to technical changes without further notice.
 © 2016 Teleflex Incorporated. All rights reserved. 94 07 08 - 00 00 01 - REV A - 07 16 PDF



Distributed by:
Teleflex Headquarters International, Ireland · Teleflex Medical Europe Ltd. · IDA Business & Technology Park
 Dublin Road · Athlone · Co Westmeath · Tel. +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73 · orders.intl@teleflex.com
United Kingdom Tel. +44 (0)14 94 53 27 61 · info.uk@teleflex.com
South Africa Tel. +27 (0)11 807 4887 · assist.africa@teleflex.com