

WatchBP[®] office ABI

Cardiovascular Screening Device


1. Dual-cuff for simultaneous both-arm blood pressure measurement.
2. Atrial fibrillation (AFIB) detector.
3. Automated ankle brachial index (ABI) assessment function.

Features:

Unique dual-cuff design for convenient, simultaneous both-arm blood pressure measurement for assessing inter-arm blood pressure differences.

With one touch the device automatically performs duplicate measurements according to guidelines while also screening for **AFIB**.

Fast, easy and reliable automated **ABI** assessment, screening for peripheral arterial disease without Doppler.

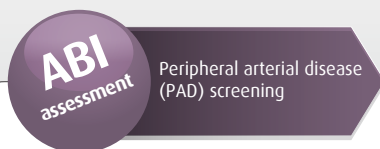

microlife[®]

Automated ABI Assessment

✓ Accurate

A clinical validation study has shown that the WatchBP Office ABI is highly accurate and that ABI values obtained with the WatchBP Office ABI are as reliable as when obtained with a Doppler device.^[4]

✓ Reliable


Sensitivity and specificity of the WatchBP Office ABI to diagnose Peripheral arterial disease (defined as Doppler ABI <0.9) at 85% and 97% respectively.^[4]

✓ Fast and Easy

Average time for automated ABI assessment is only 4 min with the WatchBP Office ABI. It requires less than half the time compared to Doppler measurement.^[4]

✓ PC Connectivity

Transmit all measurement data to any PC via USB connection.

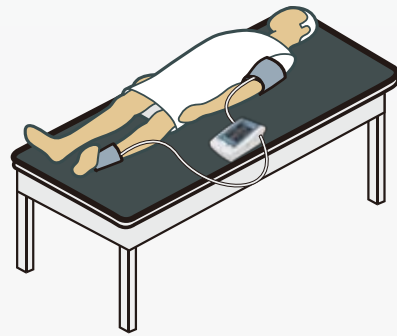

Automated AFIB Detection

✓ Accurate Screening

- Reliable detection of AFIB during automated blood pressure measurement with 97-100% sensitivity and 89% specificity as compared to ECG diagnosis.^[5,6]
- Easy, convenient and screening for AFIB function which follows the stroke prevention guidelines of leading medical society.^[7]

References

1. World Health Organization 2007. www.ish-world.com/default.aspx?Guidelines.
2. O'Brien, et al. J Hypertens 2003; 21:821-48.
3. Pickering, et al. Hypertension 2005; 45:142-61.
4. Kollias, et al. J Clin Hypertens 2010;12(Suppl 1):A81.
5. Wiesel J, et al. Am J Hypertens 2009;22:848-52.
6. Stergiou, et al. J Hum Hypertens 2009;23:654-8.
7. National Stroke Association Guidelines. www.stroke.org/site/PageServer?pagenam=PREVENT.